


Central Sterile Services Department

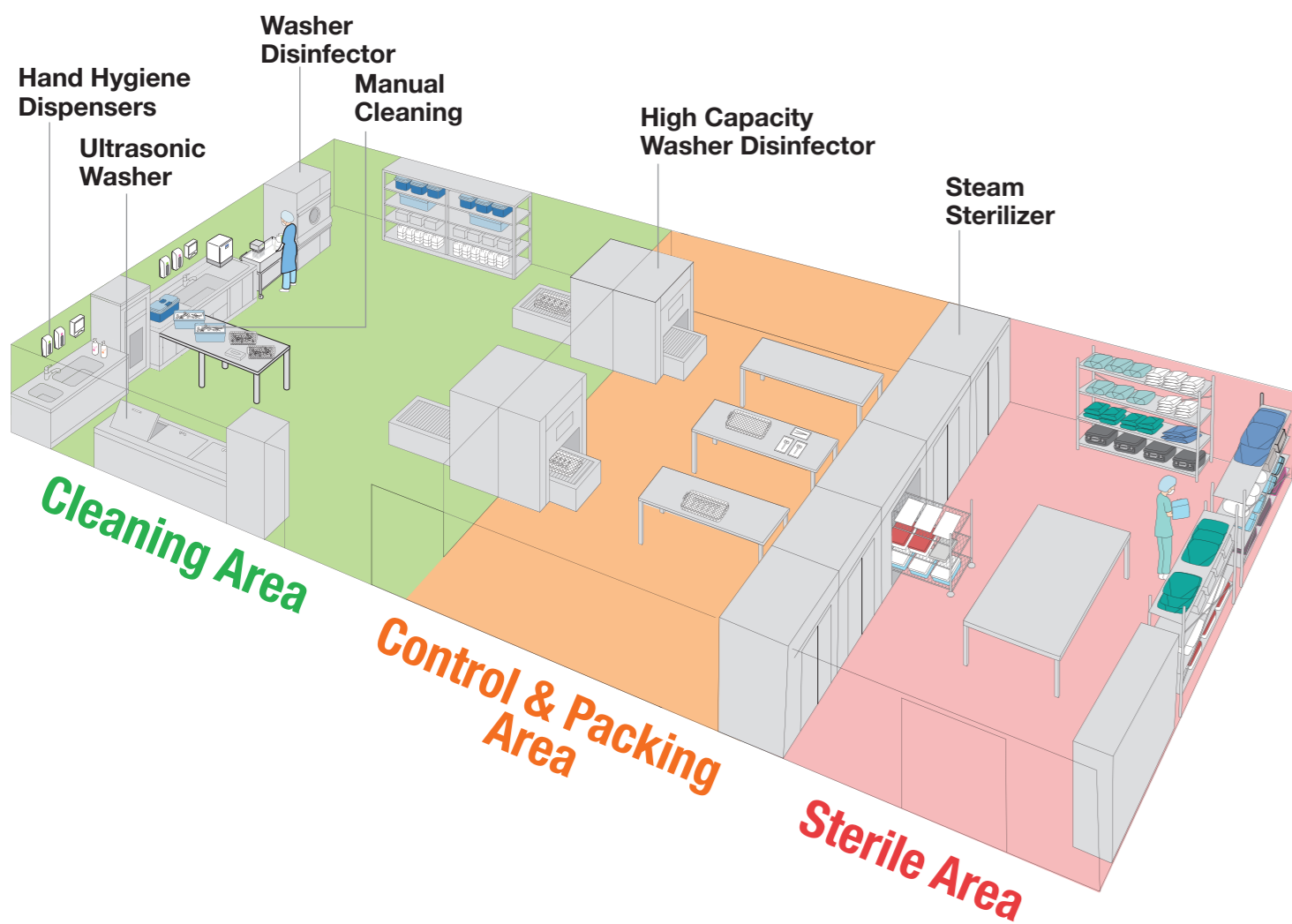


Pre-Cleaning	Pre and Main Cleaning				Main Cleaning	
Pre-Cleaning Spray	Enzyme Cleaner for Manual Soaking	Enzyme Cleaner for Manual Soaking	Multitype Enzyme Cleaner for Manual Soaking	Multi-Purpose Enzyme Cleaner	Cleaner for Automated Washers	Cleaner for Automated Jet Washers
Maintenance						
Cleaner for Ultrasonic Washers	Lubricant, Rust Inhibitor, and Drying Aid Power Forge	Heat Stain Remover	Rust Remover S	Rust Remover	Residual Protein Detector	

Medical Device Reprocessing System


	Enzyme Cleaner for Manual Soaking	Enzyme Cleaner for Manual Soaking	Multitype Enzyme Cleaner for Manual Soaking	Multi-Purpose Enzyme Cleaner	Cleaner for Automated Washers	Cleaner for Automated Jet Washers	Cleaner for Ultrasonic Washers
Cleaning method	Soaking	Soaking Ultrasonic cleaning	Soaking Ultrasonic cleaning	Soaking, Jet cleaning, Ultrasonic cleaning	Jet cleaning Ultrasonic cleaning	Jet cleaning	Ultrasonic cleaning
pH	Neutral	Neutral	Weak Alkaline	Weak Alkaline	Neutral	Alkaline	Alkaline
Dilution rate	1.0% (100mL per 10L of warm water)	1.0% (100mL per 10L of warm water)	1.0% (100mL per 10L of warm water)	● Soaking 1.0% (100mL per 10L of warm water) ● Jet cleaning 0.5% (50mL per 10L of water) ● Ultrasonic cleaning 1.0% (100mL per 10L of water)	● Jet cleaning 0.5% (50mL per 10L of water) ● Ultrasonic cleaning 1.0% (100mL per 10L of water)	0.5% (50mL per 10L of water)	0.5% (50mL per 10L of water)
Foaming	Foaming	Low foaming	Low foaming	Low foaming	Low foaming	Non-foaming	Non-foaming
Rinse	Necessary	Necessary	Necessary	Necessary	Necessary	Necessary	Necessary
Compatibility	Stainless steel, copper, brass and aluminum	Stainless steel, copper, brass and aluminum	Stainless steel, copper, brass and aluminum	Stainless steel, copper, brass and aluminum	Stainless steel, copper and brass	Stainless steel and aluminum	Stainless steel
Ideal usage	For manual cleaning	For manual cleaning and ultrasonic cleaning	For heavily soiled instruments	For jet cleaning, ultrasonic cleaning and manual cleaning	For facilities with pH regulated water	For jet cleaning	For ultrasonic cleaning

	Pre-Cleaning Spray	Lubricant, Rust Inhibitor, and Drying Aid Power Forge	Heat Stain Remover	Rust Remover S	Rust Remover	Residual Protein Detector
Cleaning method	Spraying	Soaking Jet cleaning	Soaking	Soaking	Soaking	Direct application
pH	Neutral	Neutral	Acidic	Neutral	Mild Acidic	Acidic
Dilution rate	No dilution	0.5% (50mL per 10L of water)	5% (200mL per 4L of water)	50% (500mL per 500mL of water)	50% (500mL per 500mL of water)	No dilution
Foaming	Non-foaming	Low foaming	Non-foaming	Low foaming	Low foaming	Non-foaming
Rinse	Unnecessary	Unnecessary	Necessary	Necessary	Necessary	Necessary
Compatibility	Stainless steel, copper, brass and aluminum	Stainless steel, copper, brass and aluminum	Stainless steel	Stainless steel	Stainless steel	Test before using
Ideal usage	When instruments cannot be cleaned immediately	For preventing rust and reducing drying time	For removal of heat stains	For removal of heavy rust	For removal of rust	For cleanliness check

Optional products



Thermostatic Chamber
TB1000P
Maximizing enzyme detergent efficacy.




500mL
Measuring
Cup
For measuring

Product	Content	Quantity per box	JAN Code
Power Quick Pre-Cleaning Spray	500mL	12	49-87696-50330-4
Power Quick Enzyme Cleaner for Manual Soaking Neutral	1L	6	49-87696-50333-5
	4L	3	49-87696-50334-2
Power Quick Enzyme Cleaner for Manual Soaking Neutral, Low foaming	1L	6	49-87696-50331-1
	4L	3	49-87696-50332-8
Power Quick Multitype Enzyme Cleaner for Manual Soaking Weak Alkaline, Low foaming	1L	6	49-87696-50335-9
	4L	3	49-87696-50336-6
Power Quick Multi-Purpose Enzyme Cleaner Weak Alkaline	4L	3	49-87696-50354-0
Power Quick Cleaner for Automated Washers Neutral	4L	3	49-87696-50337-3
Power Quick Cleaner for Automated Jet Washers Alkaline	4L	3	49-87696-50338-0
Power Quick Cleaner for Ultrasonic Washers Alkaline	4L	3	49-87696-50339-7
Power Quick Lubricant, Rust Inhibitor, and Drying Aid Power Forge	4L	3	49-87696-50340-3
Power Quick Heat Stain Remover	1L	6	49-87696-50341-0
Power Quick Rust Remover S	1L	6	49-87696-50355-7
Power Quick Rust Remover	4L	3	49-87696-50342-7
Power Quick Residual Protein Detector	250mL	12	49-87696-50343-4

■ For product development purposes, products are subject to change without notice.
 ■ Because print and photo qualities vary, products may appear different than depicted in this brochure.
 ■ The information in this catalog is current as of June, 2016.

PQ POWER QUICK Medical Device Reprocessing System

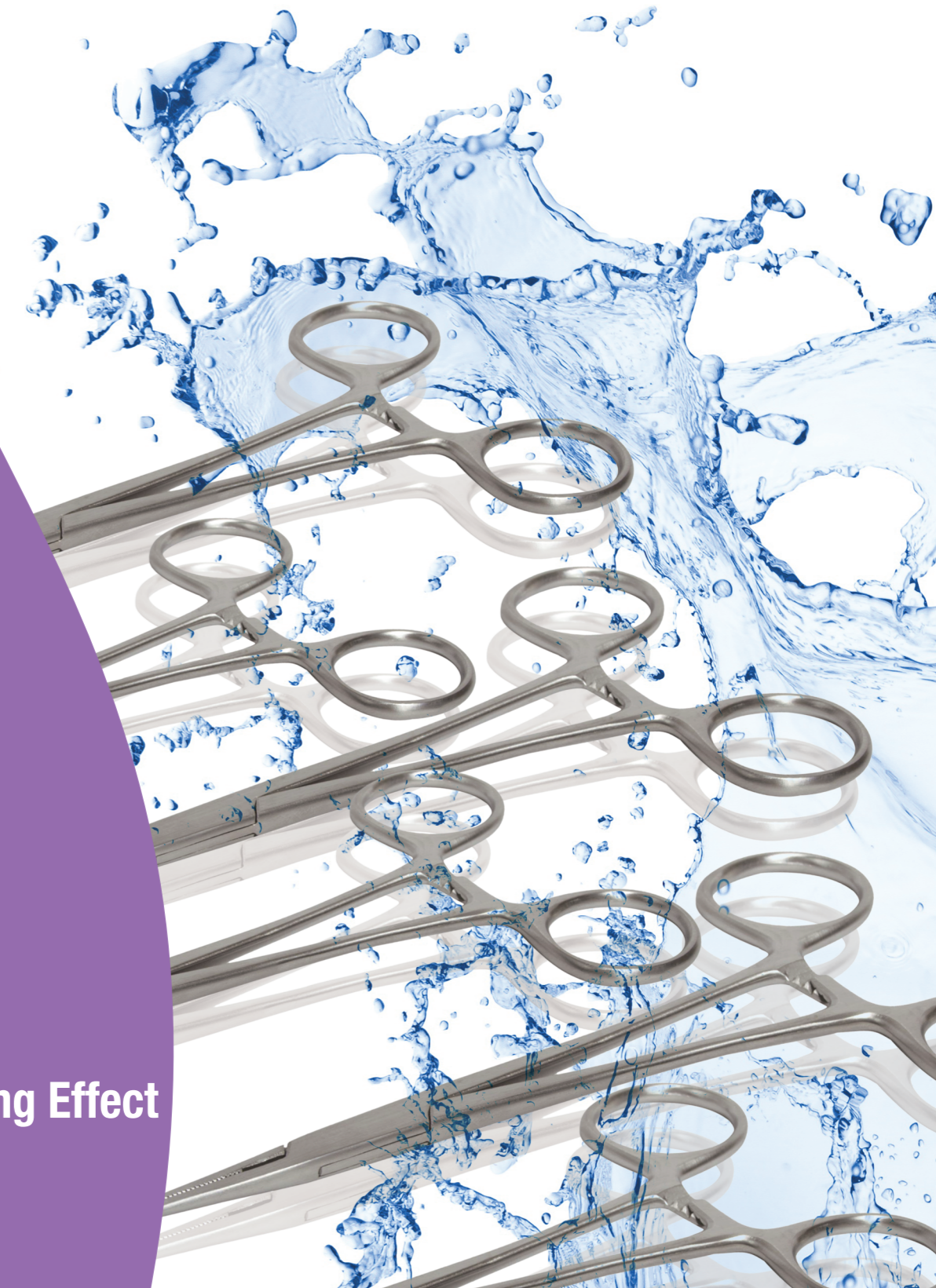

**Safety,
Efficiency,
Proven Cleaning Effect**